

# specifications

The Floor Pedestal shall be an Amico Alert-1 series.

The floor pedestal shall be made of stainless steel with a #4 satin finish.

It shall be capable of holding medical gas outlets and electrical devices within the structure of the unit.



The electrical devices are provided with electrical boxes and 18" conduit and wire.

There are (2) standard sizes to choose from:

Small Shell - 12" x 10" x 20" (305mm x 254mm x 508mm)

Large Shell - 18" x 10" x 20" (457mm x 254mm x 508mm)

The Floor Pedestal is UL Listed to U.S. and Canadian safety standards.

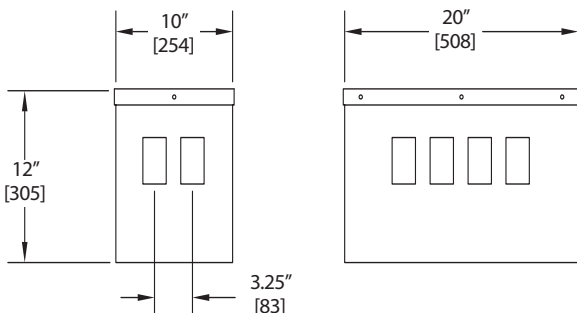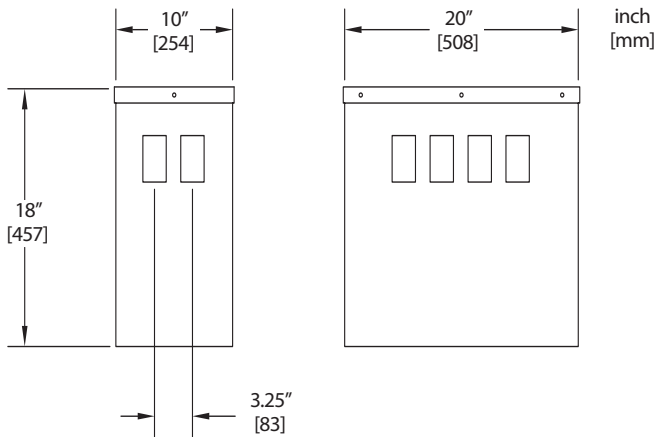

# features

- Stainless steel shell for aesthetic appearance
- All devices and gas outlets are accessible through the access panel on the back of the pedestal or the removable top
- Custom sizes are available upon request

project

---




---



---


# dimensions

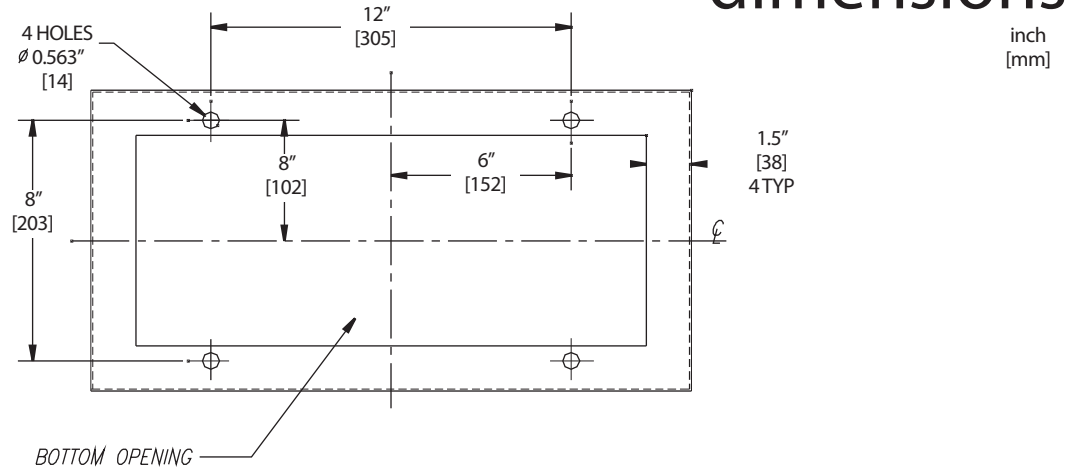


# installation

Use (four) 3/8"-1/2" diameter mounting bolts (provided by contractor) to mount the pedestal securely on the finished floor. Remove the top cover and access panel for piping and wiring.

## Small Pedestal

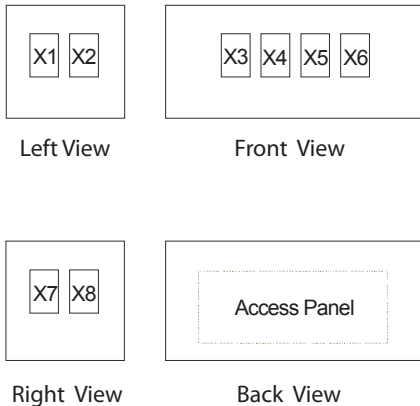

## Large Pedestal

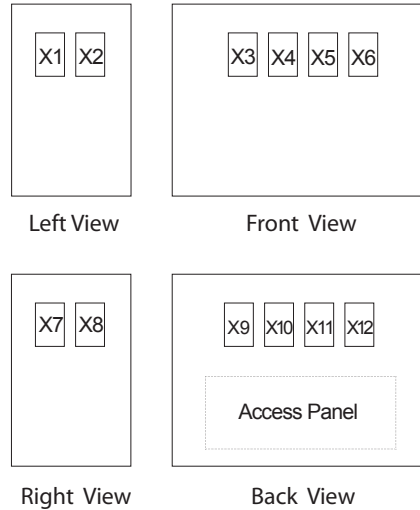

Equipment Layout	
X1	
X2	
X3	
X4	
X5	
X6	
X7	
X8	
X9	
X10	
X11	
X12	

# model numbers

Small Pedestal: W-PEDESTAL-12  
 Large Pedestal: W-PEDESTAL-18

# outlet styles

- DISS
- Puritan-Bennett
- Chemetron
- Ohmeda
- Oxequip

represented by: